FORTGESCHRITTENEN PRAKTIKUM II

Mößbauer effect

04.04.2016

Benjamin Winkelmann Peter Spalthoff

Tutor: Veronika Magerl

Contents

List of Figures							
1	Goal of the experiment						
2	physical principles						
	2.1	Interaction of Gamma radiation with matter	1				
	2.2	Gamma Decay and resonance absorption	2				
	2.3	Doppler shift	3				
	2.4	Mößbauer effect	3				
	2.5	Debye-Waller factor	4				
	2.6	Isomer shift	4				
	2.7	Hyperfine splitting	4				
3	Ехр	Experimental setup and procedure					
	3.1	Method	5				
	3.2	Setup	5				
	3.3	The source Co-57	6				
	3.4	Americium sample	6				
	3.5	Procedure	8				
4	Analysis						
	4.1	Identifying the Fe-57 transition peak	9				
	4.2		10				
	4.3		12				
	4.4		12				
5	Refe	erences	15				

List of Figures

2.1	Compton effect: A photon is scattered by a (quasi) free electron changing
	its direction by an angle $\varphi[2]$
2.2	spontaneous γ emission
2.3	Lorentz distribution
2.4	Debey factors
2.5	Hyperfine splitting Fe-57
3.1	Setup overview
3.2	Co-57 decay
3.3	Americium sample
3.4	K_{α} and K_{β} lines of Cu, Rb, Mo, Ag, Ba, and Tb[5]
3.5	The 5 visible lines of the americium source[7]
4.1	Reference spectra
4.2	MCA calibration
4.3	Compton background data
4.4	Stainless steel spectrum

1 Goal of the experiment

By measuring absorption of photons of the 14.4kev transition of Fe-57, in stainless steel and natural iron, the isomeric shift, the effective absorber thickness, the Debeye-Waller-factor of the source the lifetime of the excited Fe-57 state, the magnetic field at the location of the nucleus and the magnetic moment of the 14.4keV state.

2 physical principles

2.1 Interaction of Gamma radiation with matter

Photons interact with matter in three major ways[2]:

Photoelectric effect

Shell electrons of atoms absorb photons and gain its energy, leaving the potential well of the atom and exiting the shell with the energy $E_e = E_{\gamma} - E_B$ with E_B being the binding energy of the electron.

Compton scattering

Compton Scattering is the elastic scattering of photons at quasi free electrons ($E_B << E_{\gamma}$) and its wavelength $\lambda = 2\pi c/\omega$ is shifted, depending on the scattering angle φ (see figure 2.1):

$$\lambda_S - \lambda_0 = \frac{2\pi\hbar}{m_e c} (1 - \cos(\varphi)) \tag{2.1}$$

Pair Production

A photon can produce a positron electron pair if it has an energy of at least $2 \cdot m_e = 1.022 MeV$ the photon is lost in the processes reducing the Intensity of the photon beam. Due to those processes the intensity of electromagnetic radiation decreases exponentially with penetration depth d:

$$I(d) = I_0 \cdot exp(-\mu d) \tag{2.2}$$

where μ is the attenuation coefficient.

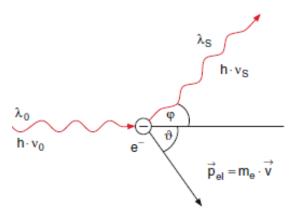


Figure 2.1: Compton effect: A photon is scattered by a (quasi) free electron changing its direction by an angle $\varphi[2]$

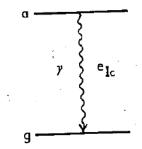


Figure 2.2: principle of spontaneous γ emission of excited nuclei. Transitioning from an excited state (E_a) to the ground state (E_g) the nucleus emits a photon with energy $E_a - E_g = \hbar \dot{\omega}$ or transmits that energy directly to an electron of the atomic shell.[1]

2.2 Gamma Decay and resonance absorption

Nuclei in excited states (energy E_a) can spontaneously transition into the ground (energy E_g) state. The energy ΔE the nucleus loses is either carried by an emitted photon (spontaneous emission) or directly gained by a shell electron (inner conversion). In the case of spontaneous emission, the photon can be absorbed by a nucleus of the same kind which thereby transits into an excited state. This is called resonance absorption. However due to the recoil the nuclei receive this rarely happens for free atoms. Consider the rest frame of a nucleus, that means its momentum is $p_0 = 0$. Now consider this

nucleus decays by emitting a photon. Since the photon carries the momentum

$$p = \frac{E_{\gamma}}{c} = \frac{\hbar}{c} \cdot \omega \tag{2.3}$$

with the reduced Planck constant \hbar and c the speed of light, and momentum is conserved the nucleus receives a recoil equal to the photons momentum and therefore also kinetic

energy or recoil energy R. The emitted photon therefore has the energy [4]:

$$E_{\gamma} = \Delta E - \frac{p^2}{2m} = \Delta E - \frac{(\hbar\omega)^2}{2mc^2} =: \Delta E - R \tag{2.4}$$

When the photon is absorbed the same applies: the absorbing nucleus receives the recoil R. This means for a photon to be absorbed, inducing a nuclear transition with ΔE the photon has to have the energy:

$$E_{\gamma} = \Delta E + R \tag{2.5}$$

In consequence the absorption spectrum is shifted relative to the emission spectrum (see fig 2.3 depending on the recoil energy R.

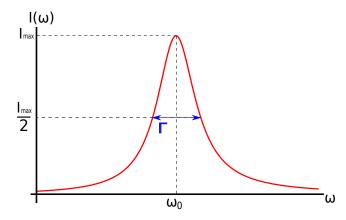


Figure 2.3: Lorentz distribution: $I(\omega) \propto \frac{1}{(\omega - \omega_0)^2 + (\Gamma/2)^2}$

2.3 Doppler shift

Due to thermal motion the emitting nucleus and the absorbing nucleus have relative velocity v, shifting the frequency via Doppler effect:

$$E_{\gamma}' = E_{\gamma}(1 + \frac{v}{c}) \tag{2.6}$$

So the energy is changed by:

$$E_{\gamma}' - E_{\gamma} = E_{\gamma} \frac{v}{c} \tag{2.7}$$

2.4 Mößbauer effect

The Mößbauer effect is the name for the phenomenon of recoilless emission (or absorption). Revisiting equation 2.4 one an see that the recoil energy $R = \frac{(\hbar\omega)^2}{2mc^2}$ is inversely proportional to the mass of the nucleus. In a solid it is possible for the whole lattice the absorb the recoil therefore increasing the recoiled mass enormously, so that $R \approx 0^1$. This effect means atoms that shows this behavior (In this experiment Fe-57) can emit photons, that can be reabsorbed by atoms of the same kind.

 $^{^{1}{\}rm this}$ is a simplified description, for a more detailed one see [4] and [1]

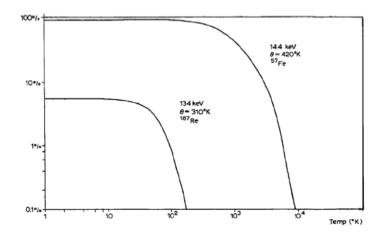


Figure 2.4: Debye factor as a function of temperature in Fe-57 and Re-187. At room temperature Fe-57 has a ratio of recoilless emission absorption of 0.91

2.5 Debye-Waller factor

The Debye-Waller factor is the ratio of recoilless absorption.

$$f = exp \left[-\frac{3R}{2k \cdot \Theta_D} \left(1 + \frac{4T^2}{\Theta_D^2} \int_0^{\Theta/T} \frac{x dx}{e^x - 1} \right) \right]$$
 (2.8)

Where R is the recoil energy, k the Boltzmann constant, Θ_D the Debye temperature. If the temperature is low $T \ll \Theta_D$ this can be simplified to:

$$f \approx exp\left[-\frac{3R}{k \cdot \Theta_D} \left(\frac{3}{2} + \frac{\pi^2 T^2}{\Theta_D^2}\right)\right]$$
 (2.9)

2.6 Isomer shift

Since electrons of an atomic shell are kept within the coulomb potential of the nucleus their potential energy depends on the charge distribution in the nucleus. Transitioning to an excited state affects this distribution therefore also affecting the potential energy of the electrons. This change in energy shifts the frequency that an absorbed photon must have to induce the transition[1].

2.7 Hyperfine splitting

The nucleus has magnetic moment (μ_I) and spin I. In a surrounding magnetic field H the energy level splits into 2I + 1 sub-energy levels. The sub-sates are characterized by the magnetic quantum number $m_I = -I, I + 1..., I - 1, I$ and the energy difference induced is:

$$E_{HFS} = \frac{\mu_I m_I H}{I} \tag{2.10}$$

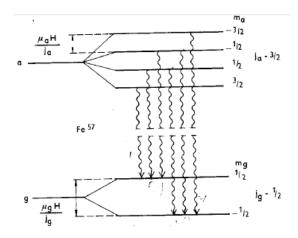


Figure 2.5: Hyperfine splitting for Fe-57 in magnetic field H for the ground state g $(I_g = 1/2)$ and the excited state $(I_a = 3/2).[1]$

The transition $m_a \to m_g$ emits photons with energy ω :

$$E_{\gamma}(m_a, m_g) = E_0 - \left(\frac{\mu_a m_a}{I_a} - \frac{\mu_g m_g}{I_g}\right) H$$
 (2.11)

3 Experimental setup and procedure

3.1 Method

To measure the absorption spectra of stainless steel and natural iron, we irradiate the samples with the 14.4keV γ -radiation emitted by a radioactive source. To vary the frequency a motor is used to move the absorber relative to the source (Doppler shift see 2.7). By repeating this measurement for different absorber velocities a spectrum is recorded.

3.2 Setup

The setup consists of the γ source, the absorber on a track, the motor used to move the absorber at constant speeds relative to the source and as the photon detector a scintillator is used. The light signal of the scintillator turned into an electric signal by a photomultiplier. This signal is amplified and shaped in the amplifier. The amplifier has two exits, one of which is connected to a single channel analyzer (SCA). If the signal pulse is within an adjustable window the SCA sends a standardized signal and enables the linear gate, which is also connected to the amplifier via a delay to ensure simultaneity of the signals. If the linear gate is enabled when it receives a signal from the amplifier it transmits the amplifier signal to the multichannel analyzer (MCA), which is read out with a Computer. The second output of the SCA is connected to a counter, which also can be read out with the Computer.

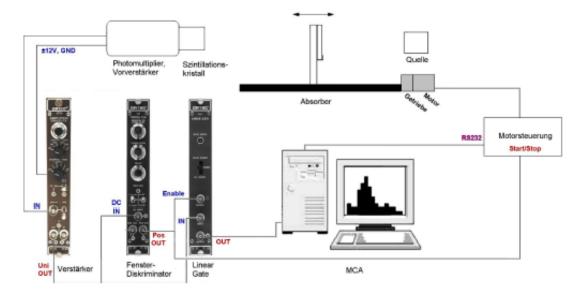


Figure 3.1: Overview of the experimental setup

3.3 The source Co-57

⁵⁷Co decays via electron capture with a branching ratio of 99.8% and a half life of 270d into an iron in an excited state ⁵⁷Fe*. This state decays with a half life of 9ns an branching ratio of 88% into the 14.4keV excited state which finally decays to the ground state (Branching ratio for $\gamma - decay$ is 10%) 3.2.

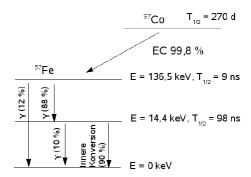


Figure 3.2: Decay series of Cobalt-57

3.4 Americium sample

To calibrate the MCA multiple reference samples are used. For this purpose americium is used as a primary source. The americium is shielded in an stainless steel case with a aperture. Attached in front of this aperture a disc with multiple targets(Cu, Rb, Mo, Ag,

Ba, and Tb) by rotating the disc one can choose a target (see fig 3.3). The radiation of the Americium source excites the target material which in turn starts emitting characteristic x-rays (x-ray fluorescence)[7]. The characteristic lines of the target samples can be found in fig 3.4 and for the americium source in fig

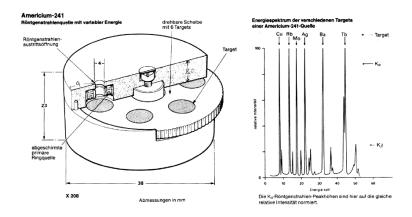


Figure 3.3: Americium sample with target revolver used as reference for the MCA calibration[5]

Target	Energi	e [keV]	Ausbeute*
	K_{α}	K_{β}	$[(\gamma/s)/sr]$
Cu	8,04	8,91	$2, 5 \cdot 10^3$
Rb	13,37	14,97	$8, 8 \cdot 10^{3}$
Mo	17,44	19,63	$2,43\cdot 10^4$
Ag	22,10	24,99	$3,85\cdot 10^4$
Ba	32,06	36,55	$4,65\cdot 10^4$
Tb	44,23	50,65	$7,6\cdot 10^4$

Figure 3.4: K_{α} and K_{β} lines of Cu, Rb, Mo, Ag, Ba, and Tb[5]

Energie [keV]	Häufigkeit [%]	Zerfallsmodus
13.927	13.0	L _o -Übergang
17.611	20.2	L_{β} -Übergang
20.997	5.2	L _γ -Übergang
26.345	2.4	E1-Kernübergang
59.536	35.7	E1-Kernübergang

Figure 3.5: The 5 visible lines of the americium source[7]

3.5 Procedure

3.5.1 MCA calibration

First the window size was set to maximum and the spectrum of the cobalt source was recorded. To identify the 14.4keV peak of the source, the (known) spectra of Cu, Rb, Mo, Ag, Ba, and Tb are measured for 300s each. The results are used to identify the 14.4keV peak in the source spectrum. The Window of the SCA was adjusted accordingly, by recording the source spectrum while adjusting the window and repeatedly resetting the recording on the computer. The Window was then adjusted until only the channels of the 14.4keV peak get a signal. We chose the settings:

upper level: 1.10lower level: 0.69

3.5.2 background measurement

The main source of background are photons of the transition between the 136.5keV state and the 14.4kev state (see fig 3.2) being scatter via Compton scattering in the scintillator and falling in the adjusted SCA window. To measure this background, aluminum plates of different thicknesses (measured with) are used to shield the scintillator. For each plate the event counts were measured over 600s. The plate thicknesses were measured with a caliber.

3.5.3 Absorption spectra of stainless steel

First a rough measurement is made: the absorption was measured for velocities of 0.1 mm/s to 1.1 mm/s (both directions) in steps of 0.1mm/s for 180s. For the finer measurements a measuring time of 300s was chosen.

3.5.4 Absorption spectra of natural iron

For natural iron the absorption was measured for absorber speeds between 0.1mm/s and 8mm/s in steps of 0.1mm/s. In a second measurement the range 0.05mm/s to 6.05mm/s was taken, also in steps of 0.1mm/s. The measuring for each velocity was 300s.

3.5.5 Attenuation through acrylic glass

The absorber is removed from the setup at the counting rate measured for 900s, once with acrylic glass and once without.

4 Analysis

4.1 Identifying the Fe-57 transition peak

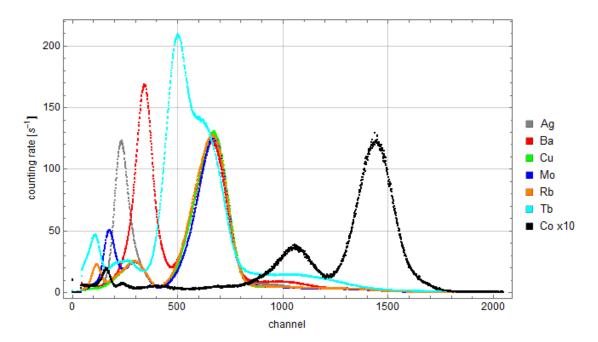


Figure 4.1: Plot of all recorded reference spectra and the cobalt source (upscaled by factor 10 for better comparability)

Around channel 700 all target samples have clearly defined peak. For Terbium(Tb) this peaks blends with its K-line. This peak is caused by photons of the americium source, passing through the targets without interacting. The K_{α} -line of copper (Cu) is beyond the left edge and therefore not measured. The peak positions are estimated from the fig 4.1. The error is estimated to be $s_{Ch} = 10$. The Result can be seen in table 4.1. One can already identify the peak since it has to lie between the Rb-peak and the Mo-peak, which suggesting the peak furthest to the left (around channel 160).

traget	energy [keV]	peak channel
Cu	8.04	-
Rb	13.37	120
Mo	17.44	180
Ag	22.10	230
Ba	32.06	340
Tb	44.23	500

Table 4.1: number of the channel for the K-line peaks

The Linear function $E(ch) = a \cdot ch + b$ was fitted to the data see figure

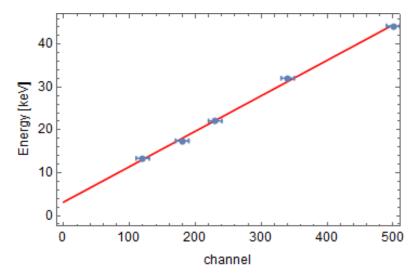


Figure 4.2: Linear fit for the calibration of the MCA

The results are:

$$a = (0.083 \pm 0.002)keV$$

$$b = (3.2 \pm 0.6)keV$$
(4.1)

The position of the 14.4keV peak should be at²:

$$ch_{14.4keV} = \frac{14.4keV - b}{a} = 135 \pm 4 \tag{4.2}$$

A Comparison with figure 4.1 shows that the peak furthest to left is the closest ($ch = 160 \pm 10$) to this value, however the values lie more almost two standard deviations apart. Fortunately, an exact calibration is not needed for further analysis.

4.2 Compton

High energy photons such as those from the transitions with 122 keV and 136 keV lose energy when passing through matter via the Compton effect. Some of these photons will randomly fall within the energy windows that was set for the measurements and thus pose an underground that needs to be subtracted from the data for parts of the further analysis. Since higher energy photons lose their energy slower than lower energy photons when passing through matter, one can separate the two by placing aluminum plates in the beam with a range of thicknesses. The count rate decreases exponentially with the thickness d

$$\dot{N} = \dot{N}_0 e^{-\mu d} \tag{4.3}$$

²error calculated according to gaussian error propagation. Unless specified otherwise, all errors of values calculated from other with other values are determined this way

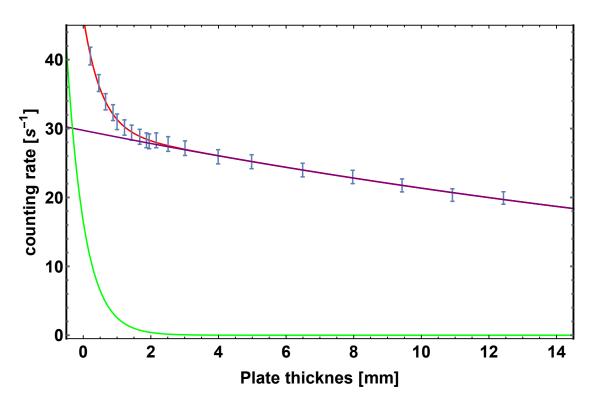


Figure 4.3: The compton background data along with the fitted sum of the two exponential functions in red. The curve for Compton underground is colored purple and that for the 14,4 keV photons green.

measurements were taken for absorber thicknesses between $d=0.21\,\mathrm{mm}$ and $d=12.43\,\mathrm{mm}$, which were measured with a caliper to such high precision that the error is far smaller than the Poisson error on the count rates and is thus neglected. Figure 4.3 shows the measured data. As two processes with different speeds, μ_C of the Compton background and μ_0 of the actual data, are expected, the fit function for the count rates for varying absorbers is

$$\dot{N}(d) = A_C \cdot e^{-\mu_C d} + A_0 \cdot e^{-\mu_0 d} \tag{4.4}$$

The resulting fit parameters for the Compton underground were

$$A_C = (29.74 \pm 0.17) \text{ s}^{-1}$$
 (4.5)

$$\mu_C = (0.0331 \pm 0.0009) \text{ mm}^{-1}$$
(4.6)

As no aluminum plates are used during regular measurements, the value for d=0, and thus the fit parameter A_C , is the underground rate to be deducted from future measurements.

4.3 Attenuation of the acrylic glass

The acrylic glass in which the sample is cased absorbs some of the radiation. To quantify this, a plate of acrylic glass of roughly the same thickness as the one used to case the sample can be but into the beam. One measurement is taken with the plate and one without. The thickness of the plate was measured as $d = (1.94 \pm 0.01)$ mm. The resulting count rates were

$$\dot{N}_0 = (107.4 \pm 0.3) \,\mathrm{s}^{-1}$$
 (4.7)

$$\dot{N}_A = (84.8 \pm 0.3) \text{ s}^{-1} \tag{4.8}$$

where the uncertainties stem from the Poisson errors on the number of counts. With the values given for the attenuation coefficient $\mu/\rho = 1.101 \text{ cm}^2/\text{g}$ and the density of the acrylic glass $\rho = 1.19 \text{ g/cm}^3$ in [5], the mass attenuation coefficient $\mu = 1.31019 \text{ 1/cm}$ can be calculated. With this value as well as the thickness of the plate, the count rate after the acrylic glass can also directly be calculated from the count rate without the acrylic glass as

$$\dot{N}_A^{calc} = \dot{N}_0 \cdot e^{-\mu d} = (83.3 \pm 0.3) \text{ s}^{-1}$$
 (4.9)

The two values agree only with their 2σ intervals. The literature value $\mu/\rho = 1.101~{\rm cm}^2/{\rm g}$ is actually listed for energies of $E_{\gamma} = 15~{\rm keV}$, which is slightly more than the that of the 14.4 keV transition. However, the value is higher for lower energies, which means that using a value for the exact energies of the photons in this experiment would lead to an even lower calculated count rate. The likely cause for the disagreement of the values is thus likely that the thicknesses of the plates do no match. A quick calculation reveals that a disagreement of 5% would bridge the gap and have the values agree within their 1σ intervals.

4.4 single line absorber

The results for the measurement of the absorption spectrum of stainless steel can be seen in figure 4.4. To evaluate the data a voigt function (convolution of a Gaussian and Lorentz) was fitted to it:

$$f(v) = B - A \cdot Voigt((v - v_0), \delta, \sigma) \tag{4.10}$$

where v is the absorber velocity, v_0 the peak position, δ is the

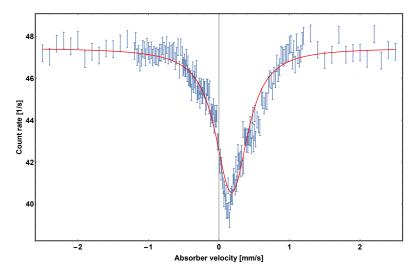


Figure 4.4: Measured absorption spectrum for stainless steel and Voigt fit.

The fit results are:

	Estimate	Standard Error
v_0	0.179	0.005
σ	0.25	0.03
δ	0.07	0.05
В	0.179 0.25 0.07 47.46	0.11
A	5.8	0.4

4.4.1 Isomeric shift

From the parameter v_0 one can calculate the isomeric shift with equation 2.7:

$$E_{iso} = (8.614 \pm 0.24) * 10^{-9} eV (4.11)$$

4.4.2 Effective absorber thickness

The effective absorber thickness is given by [5]:

$$T_A = f_A n_A \beta \sigma_0 d_A \tag{4.12}$$

where $d_A = 25\mu m$ is the absorber thickness, n_A is the number of iron atoms per volume, $\beta = 0.022$ the ratio of Fe-57 in the isotope mixture, σ_0 the cross section and $f_A = 0.8$ the Debye-Waller-factor of the absorber. To calculate n_A literature values are used ([8]):

Mößbauer effect 4 Analysis 14

$$N_A=6.022\cdot 10^{23} mol^{-1}$$
 Avogadro constant $ho_{Fe}=7.874 g/cm^3$ density of iron $A_{Fe}=55.845 g/mol$ molar mass of iron $r=(0.70\pm 0.05)$ fraction of iron in the absorber[5]

The number of Fe atoms per volume is then given by:

$$\frac{\rho_{Fe}}{A_F e} N_A r = (5.9 \pm 0.4) \cdot 10^{22} cm^{-3}$$

The absorption cross section is given by [1]:

$$\sigma_0 = \frac{\lambda^2}{2\pi} \left(\frac{2I^* + 1}{2I + 1} \right) \frac{1}{1 + \alpha} \tag{4.13}$$

where $\alpha = 9$ ([1]) is the conversion coefficient, λ the wavelength of the absorbed photon, $I^* = 3/2$ and I = 1/2 are the nuclear spins of the excited and ground state. With $\lambda = \frac{hc}{E_{\gamma}} = 86.14 \cdot pm$, the cross section is:

$$\sigma_0 = 2.36 \cdot *10^- 22m^2 \tag{4.14}$$

Plugging those values in equation 4.12, the effective absorber thickness is:

5 References

- [1] Wegener, Horst. "Der Mößbauer Effekt und seine Anwendungen". Mannheim 1966
- [2] Demtröder, Wolfgang. Experimentalphysik 3 Atome, Moleküle und Festkörper
- [3] Jakobs, Karl. Experimentelle Methoden der Teilchenphysik. Vorlesungsskript 2014
- [4] Eyges, Leonard. Physics of the Mössbauer effect. 1965
- [5] A.Zwerger (2007), S.Winkelmann (1/2011), M.Köhli (2/2011). Versuchsanleitung Fortgeschrittenen Praktikum Teil II Der Mößbauer-Effekt. 2012
- [6] M.Köhler(8/2010), M.köhli (4/2011). Versuchanleitung Fortgeschrittenen Praktikum Teil I - Kurze Halbwertszeiten
- [7] U.Landgraf. Eichspektrum der Americium-Quelle (1997)
- [8] https://www.webelements.com/ (28.4.2016)
- [9] S. Margulies, P. Debrunner, H. Frauenfelder. Transmission and line broadening in the Mössbauer effect II. 1962